

London Cancer

Central Line Care Guidelines

Based Upon UCLH Central Venous Catheter Care
Guidelines (Cancer Services) March 2010

Written by Liz Simcock

Updated and reviewed by Simon Jenkinson and Liz Simcock
for the Chemotherapy ERG May 2014

Review date: May 2016

Contents

Preface and Acknowledgements / Links	3
 Chapter 1: Care of Non-Tunnelled CVCs	4
 Chapter 2: Care of Tunnelled CVCs (Hickman Lines)	7
 Chapter 3: Care of PICCs	10
 Chapter 4: Care of Implantable Ports (Portacaths)	13
 Chapter 5: Care of CVCs used for Blood Processing (eg Apheresis)	17
 Chapter 6: Managing Complications	18
 Chapter 7: Maintaining Patency	21
 Appendix I: Background Information: Different Types of CVC	24
 Appendix II: Background Information: Principles of Care	30
 Appendix III: Background Information: Complications of CVCs	35
References	40

Preface and Acknowledgements

These guidelines are designed as a resource for nurses caring for patients with various forms of Central Venous Catheter (CVC). They were created at UCL Hospitals NHS Foundation Trust and have been adopted as a London Cancer generic guideline. They are intended to deal primarily with practical matters. No attempt has been made to address the many other issues which may face patients with CVCs - e.g. anxiety or altered body image.

Many of the procedures referred to in this guideline can be found at:

www.uclh.nhs.uk/cvc

This version of the guidelines includes the new *Tegaderm*TM dressings introduced in July 2010.

Links

These guidelines are part of a package of resources designed to help nurses caring for patients. Other resources include:

UCLH external website

www.uclh.nhs.uk/cvc

Guidelines on many of the described procedures are found on this website.

Other important links

Department of Health Saving Lives Campaign: www.clean-safe-care.nhs.uk

Chapter 1: Care of Non-Tunnelled CVCs

Infection Prevention & Management	Non-Tunnelled CVCs
	<ul style="list-style-type: none">• Infection is the most common complication associated with CVCs and one of the most serious.• You must be familiar with the <u>infection control guidance</u> when caring for patients with CVCs refer to your local hospital policy.• Decontaminate hands before and after each patient contact using correct <u>hand hygiene</u> procedure.• Always use ANTT when accessing the CVC.• Regularly inspect for signs of infection (at least daily if patient is in hospital)• Take action immediately if there are signs of CVC-related infection. These include:<ul style="list-style-type: none">○ Pyrexia, rigor, malaise (see Chapter 6 Managing Complications)○ Tenderness, inflammation and or pain at exit site (see Chapter 6 Managing Complications)• Central Venous Catheters should be removed as soon as possible if they are not needed.
Assessing Patency	Non-Tunnelled CVCs
	<p>NB Always use ANTT when accessing the CVC</p> <ul style="list-style-type: none">• Do not administer chemotherapy, drugs or fluids unless the line is fully patent. By <i>fully patent</i> we mean that:<ul style="list-style-type: none">○ The line can be flushed easily○ There is flashback of blood• If the catheter is not fully patent see Chapter 7 Maintaining Patency.• Testing for patency:<ul style="list-style-type: none">○ Test for flashback of blood before administering IV medication but note that you should not discard blood unnecessarily. To assess for flashback you can either:<ul style="list-style-type: none">▪ attach a syringe containing 10mls 0.9% saline* to the catheter, flush a couple of mls into the line and then withdraw. As soon as you see a trace of blood in the catheter or syringe just flush the rest of the saline* into the line.▪ or use a gravity technique (ie with all clamps open briefly hold an attached infusion below the level of the patient's heart until you see flashback of blood).○ BUT NB: if there are infusional vasoactive drugs in the lumen, withdraw prior to flushing to avoid bolus dose.
Flushing	Non-Tunnelled CVCs
	<p>NB Always use ANTT when accessing the CVC</p> <ul style="list-style-type: none">• Before flushing<ul style="list-style-type: none">○ If there are infusional vasoactive drugs in lumen, withdraw prior to flushing to avoid bolus dose.• Technique:<ul style="list-style-type: none">○ Brisk push-pause technique with positive pressure finish• What to flush with:<ul style="list-style-type: none">○ 10 mls 0.9% saline* between incompatible drugs / infusions and after blood sampling○ Lock with a further 10mls 0.9% saline*• Frequency of flushing:<ul style="list-style-type: none">○ Flush unused lumens at least once a week (2 x 10mls 0.9% saline*)

*must be prescribed

NB Always use ANTT for exit site care

- **Securement:**

- Always fix catheter firmly to patient's skin with *Statlock*[®] dressing.

- **Cleaning:**

- Clean exit site at dressing changes with *Chloraprep 3mIs*[®] using a 30-second back and forth friction rub. Allow to dry.
- NB If there is loose blood or exudate present, this should be removed first using sterile gauze and 0.9% sterile saline.



- **Dressings:**

- Patients aged 16+: *Statlock*[®] dressing plus *Tegaderm*[™] CHG or equivalent dressing (right). Change every 7 days (or sooner if dressing becomes wet, soiled or detached).
- Under 16s: *Statlock*[®] dressing plus a transparent dressing. Change every 7 days (or sooner if dressing becomes wet, soiled or detached).
- If patient cannot tolerate *Tegaderm*[™] CHG, try a plain transparent dressing with *Biopatch*[®]. Change every 7 days (or sooner if dressing becomes wet, soiled or detached).
- If patient cannot tolerate a transparent dressing at all, use gauze-type dressing (eg *Mepore*[®]) with *Biopatch*[®]. Inspect and change daily if patient at high risk of exit site infection. (Change sooner if dressing becomes wet, soiled or detached).



- **Bathing & showering**

- The exit site must not be allowed to get wet.

- **Central Venous Catheters** should be removed as soon as possible if they are not needed.

- **Who can remove non-tunnelled CVCs?**

- Any qualified nurse who is IV trained and has been assessed as competent and who follows these guidelines.

- **Procedure:**

- **You will need assistance** during this procedure: do not attempt it alone.
- **Check patient's coagulation status.** If there is an increased risk of bleeding discuss with medical team before proceeding. If platelets are < 50, platelets should be administered immediately prior to the procedure. If the patient is anticoagulated, this should be managed as for surgery.
- **The risk of air embolism increases if patient is dehydrated, is unable to lie flat, or has an uncontrolled cough.** Assess for these risks. Only proceed if satisfied that it is safe to do so.
- **Use ANTT** throughout.
- **Unless contraindicated (eg head injury or respiratory difficulties), lie the patient flat and tip the head of the bed downward** to reduce the risk of air embolism (except femoral catheters).
- **Remove the dressing.** If there is any sign of infection, take a swab of the exit site.
- **Remove *Statlock*[®] dressing.**
- **Ask patient to perform Valsalva's manoeuvre** (ie take a deep breath, hold it, and bear down). If patient unable to do this, remove the catheter during expiration and NEVER when the patient is breathing in, as this will increase the risk of air being sucked into the venous system.

- **Gently and swiftly pull out the catheter** and immediately apply pressure to the site using sterile gauze. The patient they can now breathe normally and the bed can be returned to the flat position.
- **Continue applying pressure to the exit site for three minutes** (or longer in cases of deranged clotting).
- **If systemic infection is suspected**, use sterile scissors to cut off the tip of the catheter and without contaminating it drop it into a dry sterile specimen pot. Send it to microbiology for culture.
- **Apply a sterile occlusive dressing** to prevent air from entering the venous system.
- **Advise the patient to stay in bed for 30 minutes to allow any bleeding to stop.**
- **During this time observe patient** for signs of haematoma (ie, swelling, pain, altered voice, airway obstruction).

Chapter 2: Care of Tunnelled CVCs (Hickman Lines)

Infection Prevention & Management

Tunnelled CVCs

- **Infection** is the most common complication associated with CVCs and one of the most serious.
- You must be familiar with the infection control guidance when caring for patients with CVCs.
- Decontaminate hands before and after each patient contact using correct hand hygiene procedure.
- Always use ANTT when accessing the CVC.
- **Regularly inspect for signs of infection** (at least daily if patient is in hospital)
- Take action immediately if there are signs of CVC-related infection. These include:
 - Pyrexia, rigor, malaise (see Chapter 6 Managing Complications)
 - Tenderness, inflammation and or pain at exit site (see Chapter 6 Managing Complications)
- **Tunnelled lines** should be removed as soon as possible if they are not needed.

Assessing Patency

Tunnelled CVCs

NB Always use ANTT when accessing the CVC

- **Do not administer chemotherapy, drugs or fluids unless the line is fully patent. By *fully patent* we mean that:**
 - The line can be flushed easily
 - There is flashback of blood
- **If the catheter is not fully patent see Chapter 7 Maintaining Patency** but note that community nurses administering a weekly flush to an unused lumen are not routinely required to check for flashback.
- **Testing for patency:**
 - Test for flashback of blood before administering IV medication but note that you should not discard blood unnecessarily. To assess for flashback you can either:
 - attach a syringe containing 10mls 0.9% saline* to the catheter, flush a couple of mls into the line and then withdraw. As soon as you see a trace of blood in the catheter or syringe just flush the rest of the saline* into the line.
 - or use a gravity technique (ie with all clamps open briefly hold an attached infusion below the level of the patient's heart until you see flashback of blood). NB this may not work for valved PICCs (eg *PASV*[®] or *Groshong*[®]).
 - BUT note that if there are infusional vasoactive drugs in the lumen, withdraw prior to flushing to avoid bolus dose.

Flushing

Tunnelled CVCs

NB Always use ANTT when accessing the CVC

- **Before flushing**
 - If there are infusional vasoactive drugs in the lumen, withdraw prior to flushing to avoid bolus dose.
- **Technique:**
 - Brisk push-pause technique with positive pressure finish
- **What to flush with:**
 - 10 mls 0.9% saline* between incompatible drugs / infusions and after blood sampling
 - Lock with a further 10mls 0.9% saline*.
- **Frequency of flushing:** Flush unused lumens once or twice a week (2 x 10mls 0.9% saline*).

*must be prescribed

NB Always use ANTT for exit site care.

• **Securement:**

- In addition to the exit site dressing, always fix catheter firmly to patient's skin (eg using tape or dedicated device eg *Skinfix*.)

• **Sutures:**

- **Exit site:** remove at 21 days
- **Venepuncture site:** Remove stitches / Steri-strips at 7 days (unless dissolvable)

• **Cleaning:**

- Clean exit site at dressing changes using *Chloraprep*[®] 3mls using a 30-second back and forth friction rub. Allow to dry
- NB If there is loose blood or exudate present, this should be removed first using sterile gauze and 0.9% sterile saline.



• **Dressings:**

○ **Exit site:**

- **Post-insertion:** gauze under transparent dressing for 1 day or until bleeding stops.

○ **After 1 day:**

- Patients aged 16+: *Tegaderm*TM CHG or equivalent dressing. (Under 16s: transparent dressing only.) Change every 7 days (or sooner if dressing becomes wet, soiled or detached).
- If patient cannot tolerate *Tegaderm*TM CHG or equivalent try a plain transparent dressing (with *Biopatch*[®] if 16+). Change every 7 days (or sooner if dressing becomes wet, soiled or detached).
- If patient cannot tolerate a transparent dressing at all, use gauze-type dressing - eg *Mepore*[®] (with *Biopatch*[®] if 16+). Change every 7 days or more frequently if at high risk of exit site infection. (Change sooner if dressing becomes wet, soiled or detached).



○ **After 21 days:** choose between

- Continuing with same dressing OR no dressing if appropriate for patient

○ **Venepuncture site:**

- Dry dressing and/or transparent dressing until sutures removed / dissolve.

• **Bathing, showering & swimming**

- **Bathing:** Patient should not submerge exit site in bathwater.
- **Showering:** If transparent dressing is intact patient can shower. If patient has dry dressing or no dressing, s/he can shower after 21 days as follows:
 - Remove dry dressing (if any) immediately before or after showering
 - Dry exit site after shower using sterile gauze and non-touch technique.
 - Clean exit site as usual & apply new dressing (if any).
- **Swimming:** Not generally advised because of the infection risk but if patient wishes to swim after 21 days it may be considered depending on other patient risk factors.

Patient Education and Discharge Planning

Tunnelled CVCs

- **If patient is discharged with catheter in situ** consider what information needs to be provided to the patient and/or carer to care for the CVC. Ensure this information is provided before discharge. E.g.
<http://www.macmillan.org.uk/Cancerinformation/Cancertreatment/Treatmenttypes/Chemotherapy/Linesports/Centrallines.aspx>

Removal

Tunnelled CVCs

- **Tunnelled CVCs** should be removed as soon as possible if they are not needed.
- **Do not remove tunnelled CVCs unless you have been specifically trained to do so.**

Chapter 3: Care of PICCs

Infection Prevention & Management	PICCs
<ul style="list-style-type: none">• Infection is the most common complication associated with CVCs and one of the most serious.• You must be familiar with the <u>infection control guidance</u> when caring for patients with CVCs.• Decontaminate hands before and after each patient contact using correct <u>hand hygiene</u> procedure.• Always use <u>ANTT</u> when accessing the CVC.• Regularly inspect for signs of infection (at least daily if patient is in hospital)• Take action immediately if there are signs of CVC-related infection. These include:<ul style="list-style-type: none">○ Pyrexia, rigor, malaise (see Chapter 6 Managing Complications)○ Tenderness, inflammation and or pain at exit site (see Chapter 6 Managing Complications)• PICCs should be removed as soon as possible if they are not needed.	
General Points	PICCs
<ul style="list-style-type: none">• Assess external length of PICC before use: if it has increased since insertion see Chapter 6 Managing Complications. See PICC Care Sheet 1: How to Measure a PICC on UCLH CVC webpage. www.uclh.nhs.uk/cvc• Take care at all times not to pull PICC out. Remember there's nothing to keep the PICC in apart from the dressing.• Avoid compression to vein containing the PICC. Do not use blood pressure cuff. Any bandage / tubular dressing must be loose.• Use volumetric pump with a filtered giving set when infusing blood products to avoid blockage.• Never use PICC for administering contrast medium as this will cause the PICC to split (unless the patient has a CT-Rated PICC).	
Assessing Patency	PICCs
<p>NB Always use ANTT when accessing the CVC</p> <ul style="list-style-type: none">• Do not administer chemotherapy, drugs or fluids unless the line is fully patent. By <i>fully patent</i> we mean that:<ul style="list-style-type: none">○ The line can be flushed easily○ There is flashback of blood• If the PICC is not fully patent see Chapter 7 Maintaining Patency but note that community nurses administering a weekly flush to an unused lumen are not routinely required to check for flashback.• Testing for patency:<ul style="list-style-type: none">○ Test for flashback of blood before administering IV medication but note that you should not discard blood unnecessarily. To assess for flashback you can either:<ul style="list-style-type: none">▪ attach a syringe containing 10mls 0.9% saline* to the catheter, flush a couple of mls into the line and then withdraw. As soon as you see a trace of blood in the catheter or syringe just flush the rest of the saline* into the line.▪ or use a gravity technique (ie with all clamps open briefly hold an attached infusion below the level of the patient's heart until you see flashback of blood). NB this may not work for valved PICCs (eg <i>PASV</i>[®] or <i>Groshong</i>[®]).○ BUT if there are infusional vasoactive drugs in the lumen, withdraw prior to flushing to avoid bolus dose.	

*must be prescribed

Flushing

PICCs

NB Always use ANTT when accessing the PICC

• Before flushing

- If there are infusional vasoactive drugs in the lumen, withdraw prior to flushing to avoid bolus dose.

• Technique:

- Brisk push-pause technique with positive pressure finish

• What to flush with:

- 10 mls 0.9% saline* between incompatible drugs / infusions or after blood sampling
- Lock with a further 10mls 0.9% saline*

• Frequency of flushing:

- Flush unused lumens at least once a week (2 x 10mls 0.9% saline*). Increase to twice weekly if there are patency problems.

Exit Site Care

PICCs

NB Always use ANTT for exit site care

- Refer to PICC Care guidance on the UCLH internet site www.uclh.nhs.uk/cvc for detailed guidelines on PICC exit site care.

• Securement:

- Always fix catheter firmly to patient's skin using steri-strips, *Statlock*[®] and a transparent dressing.

• Cleaning:

- Clean exit site dressing changes with *Chloraprep Sepp*[®] (see right) using a 30-second back and forth friction rub. Allow to dry.
- There is no need to clean the actual PICC itself. This is unnecessary and risks dislodging the PICC.
- NB If there is loose blood or exudate present, this should be removed first using sterile gauze and 0.9% sterile saline.



• Dressings:

- **Post-insertion:** gauze, steristrips and *Statlock*[®] dressing under a transparent dressing for 1 day (may be longer for outpatients).
- **After 1 day:**
 - Patients aged 16+: *Statlock*[®] dressing plus *Tegaderm*[™] CHG or equivalent dressing. Change every 7 days (or sooner if dressing becomes wet, soiled or detached).
 - Under 16s: *Statlock*[®] dressing plus a transparent dressing. Change every 7 days (or sooner if dressing becomes wet, soiled or detached).
 - If patient cannot tolerate *Tegaderm*[™] CHG, try a plain transparent dressing with *Biopatch*[®]. Change every 7 days (or sooner if dressing becomes wet, soiled or detached).
 - If patient cannot tolerate a transparent dressing at all, use gauze-type dressing - eg



*must be prescribed

Mepore[®] (with *Biopatch*[®] if aged 16+). Change every 7 days or more frequently if at high risk of exit site infection. (Change sooner if dressing becomes wet, soiled or detached).

- **Bathing, showering & swimming:**

- **Bathing & Showering:** Patient should not get the dressing wet. If possible provide a waterproof covering for bathing and showering (eg *Limbo*[®] or *Bathguard*[™]).
- **Swimming:** not advised unless using completely waterproof cover (available from the central venous access team).

Patient Education & Discharge Planning

PICCs

- **If patient is discharged with catheter in situ** consider what information needs to be provided to the patient and/or carer to care for the CVC. Ensure this information is provided before discharge. E.g. <http://www.macmillan.org.uk/Cancerinformation/Cancertreatment/Treatmenttypes/Chemotherapy/Linesports/PICCLine.aspx>

Removal

PICCs

- **PICCs** should be removed as soon as possible if they are not needed.
- **Who can remove PICCs?**
 - Any qualified nurse who follows these guidelines.
- **Procedure:**
 - **Always use ANTT when removing a PICC**
 - **Patient should be sitting/lying** with the PICC exit site below the level of the heart (this will help prevent air embolism)
 - **Remove the dressing & any stitches.** (Take swab if signs of infection)
 - **Pull PICC out slowly and gently** an inch or two at a time. As each inch goes by, change the position of your hand so that your fingers are close to the exit site. This will reduce the likelihood of the catheter breaking.
 - **If you meet resistance, STOP.** Resistance may be due to venospasm. If this happens, apply warm packs to the patient's arm for about 5 minutes before resuming. If there is still resistance, call the central venous access team for advice.
 - **Once PICC is out, apply pressure to exit site** with sterile gauze for 3 minutes.
 - **If systemic infection is suspected,** use sterile scissors to cut off the tip of the catheter and without contaminating it drop it into a dry sterile specimen pot. Send it to microbiology for culture.
 - **Apply sterile occlusive dressing** to prevent air from entering the venous system.
 - **Keep wound dry for 1 to 2 days** or until healed

Chapter 4: Care of Implantable Ports (Portacaths)

Infection Prevention and Management	Ports
<ul style="list-style-type: none">• Infection is the most common complication associated with CVCs and one of the most serious.• You must be familiar with <u>infection control guidance</u> when caring for patients with CVCs.• Decontaminate hands before and after each patient contact using correct <u>hand hygiene</u> procedure.• Always use <u>ANTT</u> when accessing the CVC.• Regularly inspect for signs of infection (at least daily if patient is in hospital)• Take action immediately if there are signs of CVC-related infection. These include:<ul style="list-style-type: none">○ Pyrexia, rigor, malaise (see Chapter 6 Managing Complications)○ Tenderness, inflammation and or pain at exit site (see Chapter 6 Managing Complications)	
General Points	Ports
<ul style="list-style-type: none">• Only access port using a dedicated non-coring needle with integral extension set with clamp / stopcock.• Following insertion of the port there may be oedema and tenderness around the site. This may make accessing port painful and more difficult than usual. Ideally port should be accessed while patient is in Interventional Radiology if it is to be used immediately afterwards.• Use volumetric pump with a filtered giving set when infusing blood products to avoid blockage• If patient undergoes MRI scan, inform scanning personnel about the port.• If patient requires defibrillation do not place paddles directly over the port.• The port should never be used for power-injection of contrast medium as this may cause the catheter to split (unless the patient has a CT-Rated Port).	
Inserting the Non-coring Needle	Ports
<ul style="list-style-type: none">• Which needle?<ul style="list-style-type: none">○ Style: For infusions, a 90° non-coring needle with extension set should be used. For boluses, blood-taking and flushing, a straight non-coring needle with extension set may be used instead if preferred.○ Gauge: A 20 or 22-gauge needle will suffice for most uses including blood administration and withdrawal.○ Length: Where a 90° needle is used, the length will depend on the amount of subcutaneous tissue between the skin surface and the port. The external part of the needle should not exert pressure on the skin but equally it should not stand too proud. Hint: a 1" needle is suitable for most adult patients. Deeper or more superficial ports will require longer or shorter needles.• Technique:<ul style="list-style-type: none">○ Always use ANTT when accessing the port○ Numb skin over the port if required using topical anaesthetic (before skin prep) or subcutaneous Lignocaine 1% (after skin prep).○ Prepare skin over the port using Chloraprep 3mIs[®] and using a 30-second back and forth friction rub. Allow to dry. Do not touch the proposed needle insertion site again except with totally sterile gloves.○ Prime needle and/or giving set with 0.9% saline*.○ Put on sterile gloves if you need to palpate the port to ensure you are confident of its position (otherwise you can wear non-sterile gloves).○ Hold port firmly (eg with thumb and two fingers) and stretch skin taut during insertion of the needle to prevent the port sliding out of the way of the needle, and to reduce the risk of the port becoming	

dislodged within the subcutaneous pocket.

- **Insert needle swiftly and firmly** until it is felt to contact the back of the port.
- **Verify correct position** by flushing with 10 mls 0.9% saline* and checking for aspiration of blood.
- **If there is any local discomfort and/ or oedema** in the tissues around or over the port this may indicate incorrect position of the needle. In this case needle should be removed (see below for technique) and a fresh attempt made. (You can use the same needle for up to 2 further attempts if it has not become contaminated or damaged.)
- **If the port flushes easily without any local discomfort/oedema but there is no flashback of blood**, this suggests that needle position is correct but that the catheter itself is not fully functional. See Chapter 7 Maintaining Patency.

Assessing Patency

Ports

NB Always use ANTT when accessing the port

- **Do not administer chemotherapy, drugs or fluids unless the line is fully patent. By *fully patent* we mean that:**
 - The line can be flushed easily
 - There is flashback of blood
- **If the port is not fully patent see Chapter 7 Maintaining Patency.**
- **Testing for patency:**
 - Test for flashback of blood before administering IV medication but note that you should not discard blood unnecessarily. To assess for flashback you can either:
 - attach a syringe containing 10mls 0.9% saline* to the giving set, flush a couple of mls into the line and then withdraw. As soon as you see a trace of blood in the catheter or syringe just flush the rest of the saline* into the line.
 - or use a gravity technique (ie with all clamps open briefly hold an attached infusion below the level of the patient's heart until you see flashback of blood).
 - BUT note that if there are infusional vasoactive drugs in the lumen, withdraw prior to flushing to avoid bolus dose.

Flushing

Ports

NB Always use ANTT when accessing the CVC

- **Non-accessed ports:**
 - **Flush at least every four weeks** with 10mls 0.9% saline* and lock with 4-5 mls Heparinised saline* 100 U/ml.
- **Accessed ports:**
 - **Technique:**
 - Brisk push-pause technique
 - **What to flush with:**
 - **10 mls 0.9% saline*** between incompatible drugs / infusions or after blood sampling
 - **If needle to be removed:** lock with 5 mls Heparinised saline 100 U/ml*
 - **If needle to remain in situ:** lock with a further 10mls 0.9% saline*

*must be prescribed

Removing the Needle

Ports

NB Always use ANTT when removing the needle

• Technique:

- **Lock port** with 5mls heparinised saline 100 U/ml*. Ideally, remove needle while injecting last ml to achieve positive pressure finish but use gauze to prevent spray. NB you will need to ask the patient or a third party to inject because you will need two hands for removing the needle. (NB If this is not possible, you can achieve a positive pressure finish by clamping the infusion set while injecting the final ml of flush and then remove needle as below.)
- **Stabilise the port** with one hand during needle withdrawal to avoid trauma to tissues. Take care to avoid a needle-stick injury.
- **Apply gentle pressure** to needle site with sterile gauze until minor bleeding has ceased. A plaster may be applied if necessary / desired.

Exit Site Care

Ports

NB Always use ANTT for exit site care

• Sutures:

- To side of port: remove at 7-10 days (unless dissolvable)
- Venepuncture site: Remove at 7 - 10 days (unless dissolvable)

• Frequency of needle change:

- If port in constant use for more than a week, change needle weekly using different puncture site.

• Dressings

○ Non-accessed ports:

- No dressing or exit site care required (except immediately following insertion of the port when wound should be kept covered until stitches removed and wound healed.)

○ Accessed ports:

- Pad needle with sterile gauze if necessary and cover with transparent iv dedicated dressing. Needle site should be visible for inspection.
- Tape tubing firmly to skin to prevent pulling on the needle.
- Inspect needle entry site at least daily.
- Advise patient to report any discomfort or swelling at the puncture site immediately.

• Bathing, showering & swimming

○ Non-accessed ports:

- Patient may bath, shower or swim freely once wound has healed.

○ Accessed ports:

- **Bathing:** Patient should not submerge exit site in bathwater.
- **Showering:** Patient may shower if needle site is completely covered with an occlusive dressing, taking care not to dislodge needle.
- **Swimming:** not advised while needle is in situ.

*must be prescribed

Patient Education

Ports

- If patient is discharged with port in situ:
 - **Ensure patient is aware of care** required. Various resources can be used to education the patient. E.g. <http://www.macmillan.org.uk/Cancerinformation/Cancertreatment/Treatmenttypes/Chemotherapy/Linesports/Implantableport.aspx>
 - **Ensure patient is aware of the importance of reporting complications** and has a contact number for this purpose
 - **Make arrangements for the port to be flushed every 4 weeks if not in use. This is usually best done in Daycare.**
 - **Note that Community Nursing Staff** are rarely trained to access ports. If community staff need training in use of the port, contact the cancer unit or venous access team.
 - **Patients may wish to learn to flush their own ports.** The central venous access team may be able to assist in teaching.

Removal

Ports

- Implantable ports are usually removed in Interventional Radiology.

Chapter 5: Care of CVCs used for Blood Processing (eg Apheresis)

Infection Prevention and Management	CVCs used for Blood Processing
<ul style="list-style-type: none"> As for tunnelled or non-tunnelled CVCs, whichever applies. 	
Information About Locking Solutions	CVCs used for Blood Processing
<p>NB Always use ANTT when accessing the CVC</p> <p>Locking solutions</p> <ul style="list-style-type: none"> If the catheter is to be used again within 12 hours, lock with 2 x 10 mls 0.9% saline* using vigorous push-pause flush and positive pressure finish. If the catheter is NOT to be used again within 12 hours, flush with 10 mls 0.9% saline* as above but then lock each lumen with Heparin 1000 U/ml*. The volume of heparin should be exactly the priming volume of the lumen in question (this is marked on each lumen). Each lumen locked in this way must be clearly labelled “Locked with heparin 1000 u/ml. Withdraw this before use. Do not flush into the patient”. Labels are available from the apheresis team or the central venous access team. When the catheter comes to be used again, always withdraw and discard the indwelling heparin prior to using the catheter. If you are unable to withdraw the heparin it may be necessary to flush the catheter without withdrawing but you must never do this without discussion with the patient’s medical team. This is because in effect you will be delivering a heparin bolus which might have clinical implications for the patient. 	
Assessing patency	CVCs used for Blood Processing
<ul style="list-style-type: none"> As for tunnelled or non-tunnelled CVCs, whichever applies but note you must always remove any indwelling heparin prior to assessing patency. If you are unable to withdraw the heparin it may be necessary to flush the catheter without withdrawing but you must never do this without discussion with the patient’s medical team. This is because in effect you will be delivering a heparin bolus which might have clinical implications for the patient. 	
Flushing	CVCs used for Blood Processing
<ul style="list-style-type: none"> As for tunnelled or non-tunnelled CVCs, whichever applies but see “Information about locking solutions” above 	
Exit Site Care	CVCs used for Blood Processing
<ul style="list-style-type: none"> As for tunnelled or non-tunnelled CVCs, whichever applies. 	
Patient Education (tunnelled lines)	CVCs used for Blood Processing
<ul style="list-style-type: none"> As for tunnelled CVCs, see Chapter 2. 	
Removal	CVCs used for Blood Processing
<ul style="list-style-type: none"> As for tunnelled or non-tunnelled CVCs, whichever applies. 	

*must be prescribed

Chapter 6: Managing Complications

Pyrexia plus or minus: rigor after flushing, sore throat, generally feeling unwell, hypotension, tachycardia, shock, exit site / tunnel infection

- **Possible cause:**

- Catheter Related Blood Stream Infection

- **Management:**

- **Refer to medical staff.** May be treatable without catheter removal depending on patient's clinical status and colonising organism. Microbiology opinion should usually be sought.
- **Take blood cultures** as per UCLH or local guidelines.
- **TPR & BP.** Frequency will depend on patient's clinical status.
- **If there are signs of exit site infection** see below.

Inflammation and tenderness at the exit site / insertion site / skin tunnel / port pocket plus or minus exudate

- **Possible cause:**

- Infection

- **Management:**

- **Take a swab**
- **Refer to medical staff.** In non-tunnelled CVCs, exit site infections involving pus should probably be removed to prevent infection spreading to the blood stream. In tunnelled CVCs and PICCs, a superficial exit site infection will often resolve with antibiotics, but infections involving the skin tunnel above the cuff or a port pocket are very difficult to treat. Microbiology opinion should usually be sought.
- **Stop using TegadermTM CHG or Biopatch[®] dressing** until symptoms resolve.
- **Increase frequency of dressing change & cleaning** depending on amount of exudate.
- **4 hourly TPR & BP** if patient in hospital
- **If patient also shows signs of systemic infection,** see also above.

The catheter is sluggish, or there is no flashback of blood, or there is complete blockage .

- **Possible causes:**

- Clotted blood in catheter
- Fibrin sheath (which may be diagnosed using fluoroscopy)
- Malpositioned catheter
- Build-up of lipids (Parenteral Nutrition)
- Drug Precipitation
- NB: In Implantable Ports needle may be incorrectly positioned: check before taking any other action.

- **Management:**

- See Chapter 7: Maintaining Patency.

Pain or visible swelling when catheter is used **or** fluid leaks from exit site when catheter is flushed.

• **Possible causes:**

- Malposition of catheter
- Internal catheter fracture
- Fibrin Sheath (which may be diagnosed using fluoroscopy)
- Separation of port and catheter (Implantable ports)
- NB: In Implantable Ports needle may be incorrectly positioned: check before taking any other action.

• **Management:**

- **Stop using the catheter.**
- **Refer to the central venous access team and / or medical staff:** a malpositioned catheter should usually be removed. Internal fracture cannot be repaired. If there is a fibrin sheath severe enough to cause leakage the catheter should be removed.
- **Chemotherapy:** follow Cytotoxic Policy if extravasation occurs.
- **If catheter is fractured or faulty complete Adverse Incident Form.** Do not retain the catheter but if possible record the make and lot number on the incident form.

Leakage from external portion of catheter when flushed.

• **Possible cause:**

- External catheter fracture

• **Management:**

- **Clamp or fold catheter** between the exit site and the leak to prevent air entry. If using a clamp (eg artery forceps) pad with gauze to avoid trauma to the catheter.
- **Catheter must be repaired or removed as soon as possible.** Some catheters can be repaired if equipment & expertise available. The advisability of repair will depend on the patient's clinical status as it carries a risk of infection. Contact the IV team or central venous access team for advice.
- **Complete Adverse Incident Form.** Do not retain the catheter but if possible record the make and lot number on the incident form.

Cuff protrudes from exit site (Tunnelled CVCs)

• **Possible cause:**

- Tissues within tunnel have failed to adhere to cuff & catheter has migrated out.

• **Management:**

- **Stop** using the catheter.
- **Tape catheter** firmly to skin at exit site
- **Refer to the central venous access team and / or medical staff** for catheter removal.

Increase in external length of a PICC

• **Possible cause:**

- PICC has migrated out

• **Management:**

- **Do NOT push the catheter back in**
- **Refer to insertion documentation which may indicate how many far PICC can safely migrate out.**
- **If PICC has come out by more than indicated as acceptable on insertion documentation, consult**

the central venous access team or medical team. Examination of the post-insertion CXR may reveal whether or not a CXR will need to be carried out to check tip position.

Pain or swelling of arm, neck or shoulder, with or without distension of neck / peripheral veins

- **Possible cause:**

- Thrombosis.

- **Management:**

- **Refer to the central venous access team and / or medical staff** for investigation of suspected thrombosis. It may or may not be possible to treat thrombosis without catheter removal.
- Thrombosis and infection often occur. If patient shows signs of infection, see guidance on managing infection above.

Palpitations / Abnormal ECG

- **Possible causes:**

- Cardiac arrhythmias related to CVC

- **Management:**

- **If patient is distressed or unwell as a result of abnormal rhythm, call for urgent medical assessment** and monitor vital signs
- **PICCs:** Pulling PICC out by 1 - 2cm may resolve the problem immediately. If possible, liaise with the central venous access team to check that this will not result in too high a tip position. If there is any reason to suspect that the palpitations may have another cause, check tip position before retracting the line.

Cardiopulmonary symptoms including any of the following: respiratory distress / failure apnoea, reduced O_2 saturation levels, tachycardia, bradycardia, hypotension, pallor, cyanosis, anxiety, chest pain, loss of consciousness

- **Possible causes:**

- Pneumothorax
- Air or catheter embolism
- Pulmonary embolism
- Cardiac tamponade / pericardial effusion

- **Management:**

- **Call for medical assistance** / PERT team / resuscitation team
- **Administer O_2**
- **Monitor vital signs**

Chapter 7: Maintaining Patency

Patency Problems

- **Patency problems are common** in Central Venous Catheters and include:
 - no flashback of blood
 - sluggish flow
 - complete blockage
- **Possible causes**
 - clotted blood in the catheter (most likely cause)
 - fibrin sheath
 - malpositioned catheter
 - drug precipitation
 - build-up of lipids (parenteral nutrition)
 - incorrectly positioned needle in an implantable port: check before taking any other action.

Preventing Patency Problems: good flushing techniques

- **Use a brisk 'push-pause' flushing technique** routinely when flushing the catheter – i.e. flush briskly, pausing briefly after approximately each ml of fluid.
- **Use a "positive pressure finish" when you lock the catheter** – i.e. clamp the line while you are flushing in the final ml. If there is no clamp remove the syringe from the *Bionector*[®] (or similar) while you are still injecting the final ml.

Managing Patency Problems

- **No flashback of blood**
 - Ask the patient to take deep breaths and try different positions. Flush briskly using 10mls saline*. If this fails use a thrombolytic (see below).
 - If lipids/drug precipitation suspected consult pharmacy advice for suitable agent to dissolve occlusion
- **Catheter flow is sluggish**
 - Ask the patient to take deep breaths and try different positions. Flush briskly with 10mls saline*. If this fails use a thrombolytic (see below).
 - If lipids/drug precipitation suspected consult pharmacy advice for suitable agent to dissolve occlusion
- **Catheter is completely blocked**
 - Use a 3-way tap technique to instil thrombolytic into catheter (see below).

*must be prescribed

What is a thrombolytic?

- **A thrombolytic** is a drug capable of breaking up a thrombus.
- **Urokinase** is the thrombolytic used for unblocking CVCs in Cancer Services: use 5000 units in 2mls per lumen.
- **A thrombolytic should always be prescribed.**
- **Heparin and Hepsal are NOT thrombolytics:** they are capable only of inhibiting thrombus formation.

When should you use a thrombolytic?

Use a thrombolytic to improve patency in the following situations:

- flashback of blood is absent
- free-flow of fluids is sluggish or intermittent
- resistance is felt when flushing
- the catheter/lumen is completely blocked

What if the thrombolytic fails to restore function?

- **If a thrombolytic used correctly (see below) fails to restore function, contact the central venous access team and / or medical team.** A chest x-ray may need to be carried out to check the position of the line. If a chest x-ray shows that the catheter is correctly placed, it may be worth investigating further using fluoroscopy which may reveal a fibrin sheath.
- **If the cause could be a build up of lipids from Parenteral Nutrition or drug precipitation** consult pharmacy advice for a suitable agent to dissolve occlusion.

How to use a thrombolytic

- **Arrange prescription.** (Caution if patient's clotting is severely deranged or if high doses of an anticoagulant are being given concurrently.)
- **Draw up the thrombolytic** as per manufacturer's instruction eg for Urokinase: reconstitute 10,000 unit vial with 2mls water for injection and dilute further to 4 mls. Use 2mls (5000 units) per lumen.
- **Instil the thrombolytic into the catheter and wait 2 hours (or preferably longer if possible).** But note that if the lumen is completely blocked do NOT force the thrombolytic into the catheter: see Using a Thrombolytic in a Completely Blocked Catheter (below).
- **Assess the catheter again.** NB It is best to flush the catheter with 0.9% saline (must be prescribed) prior to attempting to obtain flashback, otherwise there is a risk of creating a further blockage in the line before you have cleared it. NB There is no need to worry that you are flushing the thrombolytic into the patient: small doses can be flushed into the patient without danger unless the patient has exceptionally deranged clotting.
- **If full function has not returned** instil the thrombolytic again and leave in for longer – several hours or overnight if possible. NB the fact that urokinase has a half life of 20 minutes once it is in the patient's system is irrelevant to how long it carries on working while still in the catheter. This explains why leaving it in longer is often more effective.
- **If the procedure fails** to restore function consider whether lipids / drug precipitation could be causing a blockage. If not, refer to the central venous access team and / or medical staff: a chest x-ray may reveal malposition of the line.

Using a Thrombolytic in a Completely Blocked Catheter.

- **Attach 3-way-tap & syringes** see right. (NB 3-way taps are now contraindicated for routine IV use but are still available for this procedure. Always use a 3-way tap *without* an extension set)

- **Open clamp** (if there is one).

- **Open stopcock** to the empty syringe and the blocked catheter.

- **Pull back on the plunger** of the empty syringe to create a vacuum in the catheter. You will need to pull quite forcibly.

- **Maintain suction** with one hand and with the other hand turn stopcock so it is closed to the empty syringe and open to the syringe containing thrombolytic, which will be sucked into the catheter. Don't worry if it seems that very little thrombolytic is sucked in: even a tiny volume will reach several cm into the catheter.

- **Leave for several hours or overnight.** DO NOT CLAMP CATHETER as this will prevent the thrombolytic from penetrating into the line.


- **After this time, assess the catheter** by attempting to flush it catheter using 0.9% saline* in a 10ml syringe. Do not use excessive force. (It is best NOT to try aspirating before flushing at this stage as you may block the catheter again).

- **If the catheter is still completely blocked**, repeat the procedure: sometimes you will need to do it several times before it works. Sometimes leaving the thrombolytic in overnight seems to help¹. Don't worry about overdosing the patient: if the catheter is blocked none of the drug will actually have been flushed into the blood stream.

- **Once the catheter can be flushed**, and only then, check for flashback. If flashback is absent, administer thrombolytic as described above.

- **If the procedure fails** despite repeated attempts consult the central venous access team and / or medical team with a view to removing the catheter.

Diagram 1.
Three-way Tap Technique

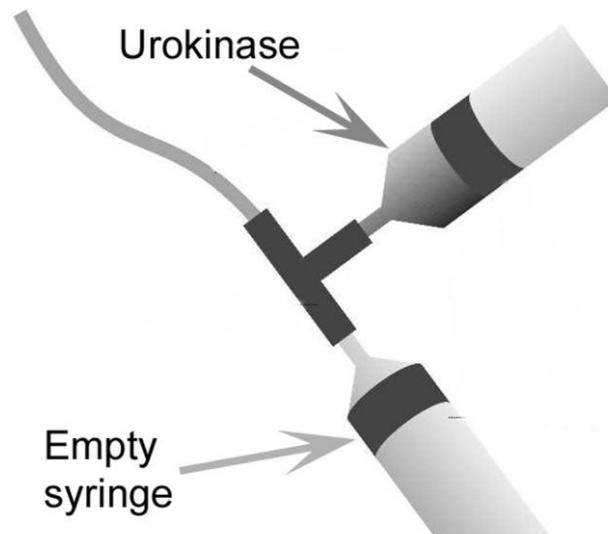

Appendix I: Background Information: Different Types of CVC

(i) Definition of a Central Venous Catheter (CVC)

The term Central Venous Catheter (CVC) refers to an **intravenous catheter whose internal tip lies in a large central vein**. There are various different types of CVC but common to all is the idea that the **tip of the catheter floats freely within the bloodstream in a large vein** and parallel to the vein wall. Blood flow around the catheter is maximised, and physical and chemical damage to the internal walls of the vein are minimised.

Opinions vary about the ideal place for the tip of a CVC² but it is generally accepted that for a catheter to be considered a “central catheter” the internal tip should be in one of the following positions.

- a. Superior vena cava
- b. Junction of the right atrium and the superior vena cava (also known as the atrio-caval junction)
- c. Right atrium
- d. Inferior vena cava above the diaphragm (femoral catheters)

Tip positions outside these areas are thought to be related to a significantly higher risk of complications, notably thrombosis^{3,4,5,6,7}.

In practice, CVC tips are not static and their position varies depending on the patient’s position, arm movements etc².


(ii) Indications

- To monitor **central venous pressure**
- To administer **large amounts** of intravenous fluids (e.g. colloids, blood products etc.)
- To administer **irritant, vesicant or hyper-osmolar** drugs / fluids (for example Noradrenaline/Adrenaline, NaHCO₃, Parenteral Nutrition, chemotherapy etc.)
- To provide **long term access** for frequent or prolonged use (e.g. chemotherapy, antibiotics, blood sampling, apheresis, continuous renal replacement therapy (CRRT), haemodialysis etc.).

(iii) Insertion and Removal

Guidelines for the insertion of central venous catheters are not covered here. **Insertion of a CVC is an invasive procedure** which must only be performed by trained, competent personnel following the latest Department of Health guidelines⁸ using “optimal aseptic technique, including a sterile gown, gloves, and a large sterile drape”⁹. Chlorhexidine 2% should be used as skinprep (except in cases of hypersensitivity). **The use of ultrasound to achieve venous access is recommended** by NICE guidelines¹⁰ but this relies upon the availability of appropriate equipment and training. Whether the catheter is inserted under general anaesthetic, sedation or simple local anaesthetic will depend upon the situation, the patient, the type of catheter to be inserted, and local practice.

Central Venous Catheters should be removed as soon as possible if they are not needed. **Techniques for the removal of a CVC** depend on the type of catheter. See Chapters 1 - 5.

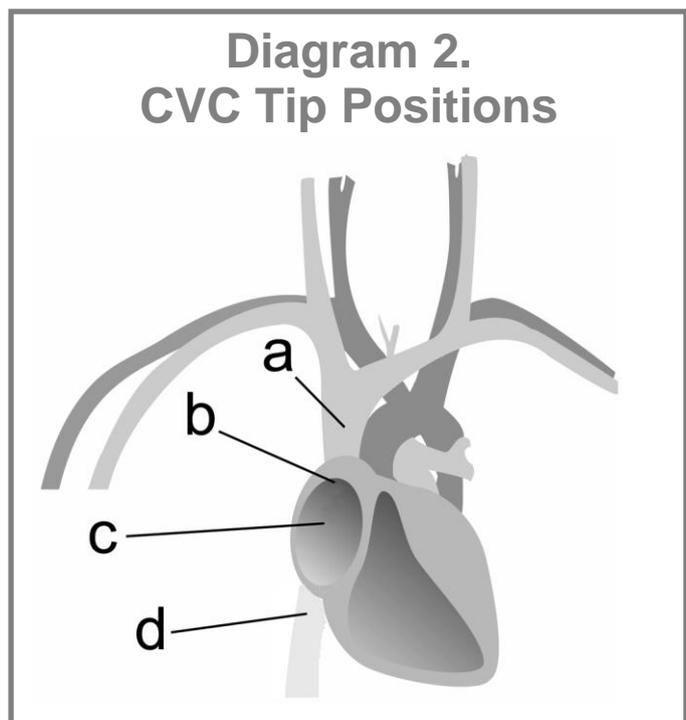


(iv) Choice of Catheter

Various different types of CVCs are available and these are described below.

The choice of device will depend chiefly on the purpose for which it is intended, though patient preference may be a key factor with long-term catheters. As a general principle **the lumen diameter and the number of lumens should be kept to a minimum**, since larger bore catheters and multiple lumens are associated with higher infection and thrombosis risks^{9,11}. Clearly there are many other factors to be weighed against these risks – e.g. in high dependency settings large bore catheters and multiple lumens tend to be used as they are essential for management of the acutely ill patient. Where parenteral nutrition is to be administered, ideally a single-lumen catheter should be used. If multiple lumens are essential, then one lumen should be dedicated “exclusively for that purpose”^{9,12}.

a) Non-Tunnelled CVCs Often called Central Lines / Neck Lines / CVP lines.


- See Chapter 1 Care of Non-Tunnelled CVCs for care of these catheters.
- Non-tunnelled CVCs are most commonly found in **acute settings**. They are not suitable for long-term use because of the risk of infection and because they are relatively uncomfortable and unsightly.
- **The catheter is usually inserted via the subclavian, jugular or femoral veins** with the tip positioned in the right atrium, the superior or inferior vena cava. It should be attached to the patient’s skin using a *Statlock*[®].
- **Non-tunnelled CVCs may have single or multiple lumens**. Each lumen provides independent access to the venous circulation, so that incompatible drugs/fluids may be administered simultaneously.
- **Each lumen is equipped with an integral clamp** to seal the catheter and guard against air entry, haemorrhage and infection.

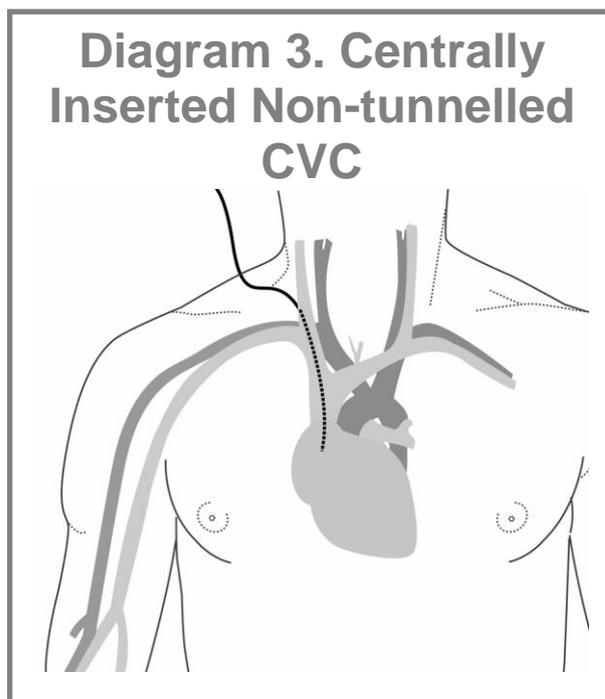

b) Tunnelled CVCs (Often called Hickman Lines)

- See Chapter 2 Care of Tunnelled CVCs (Hickman Lines) for care of these catheters.
- **Tunnelled CVCs (Hickman lines) are intended for longer-term use** in patients who require multiple infusions of fluids, blood products, drugs or parenteral nutrition. They also provide easy access for routine blood sampling. They are more comfortable and discreet than the non-tunnelled CVCs described in a) above, and can last for much longer because of the separation of the exit site and the vein (see below).
- The **tunnelled CVC** is inserted via the subclavian, jugular or femoral veins. **The catheter is tunnelled subcutaneously and exits at a convenient site** (usually on the chest wall) where it is secured with

temporary sutures. There is a 'cuff' within the tunnel to which fibrous tissue will adhere in the 2 to 3 weeks following insertion. At the end of this time, the sutures can be removed. The embedded 'cuff' helps to prevent accidental dislodgement and acts as a mechanical barrier to ascending bacteria^{13,14}.

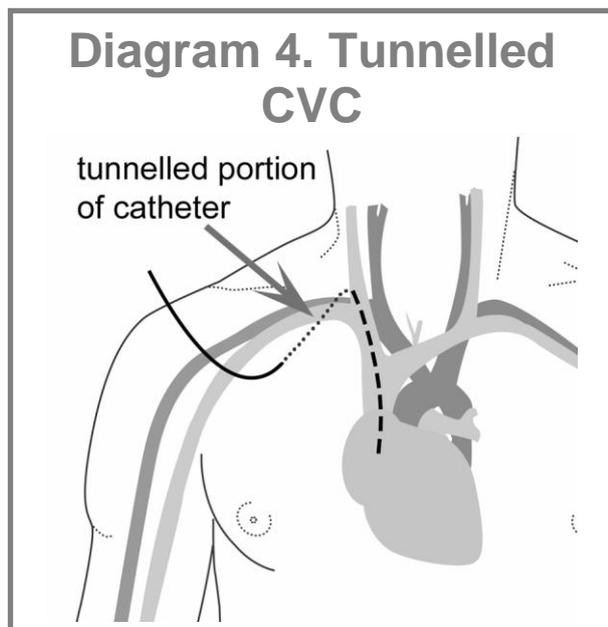


- **Single, double and triple lumen catheters are available.** Each lumen provides independent access to the venous circulation, so that incompatible drugs/fluids may be administered simultaneously.

- **Each lumen of the catheter is equipped either with an integral clamp, or a 3-way valve.** The clamp (or valve) serves to seal the catheter and guard against air entry, haemorrhage and infection. Valved catheters vary in design: the valve may be at the **internal** or **external** end of each lumen (eg *Groshong* catheters have a valve at the internal end, whereas *PASV* catheters contain a valve at the external end).

- **Patients with tunnelled CVCs may be discharged home with the catheter in situ.** In these cases patient education regarding the recognition and reporting of complications is of great importance.

Where possible, care in hospital should be aimed at the promotion of independence in caring for the line, but liaison with the primary health-care team remains vital.

c) PICCs (Peripherally Inserted Central Catheters).

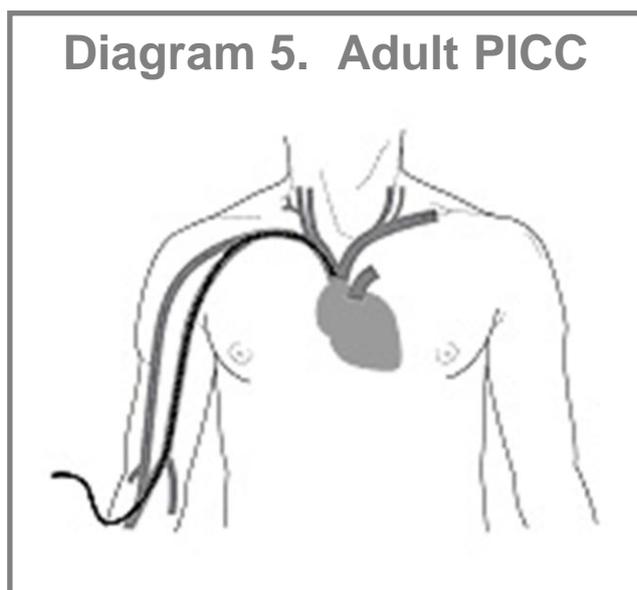

- See Chapter 3 Care of PICCs and PICC Care Page for care of these catheters.

- PICCs (Peripherally Inserted Central Catheters), like tunnelled lines, are intended for **mid to long-term use** in patients who require multiple infusions of fluids, blood products, drugs or Parenteral Nutrition.⁵⁵ They may also provide access for routine blood sampling.

- **A PICC is a fine bore CVC** inserted into a peripheral vein – usually the basilic or brachial vein – and threaded upwards towards the heart. Tip position is verified by chest x-ray following insertion (unless the tip has been screened during insertion using Fluoroscopy).

- Unlike tunnelled Lines, **PICCs do not possess a "cuff" to secure the catheter.** There is nothing to keep the PICC in place unless it is secured to the skin of the patient's arm using steri-strips and a *Statlock*[®]. In some Trusts PICCs may be stitched in position. Checking the external length of the PICC should be a routine part of care before administering drugs or fluids. See Chapter 3 and PICC care (www.uclh.nhs.uk/cvc) for detailed guidelines.


- **PICCs can be single or double lumen.** Each lumen provides independent access to the venous circulation, so that incompatible drugs/fluids may be administered simultaneously.

- **Each lumen of a PICC is equipped either with an integral clamp, or a 3-way valve.** The clamp (or valve) serves to seal the catheter and guard against air entry, haemorrhage and infection. Valved PICCs vary in design: the valve may be at the **internal** tip of each lumen (eg the *Groshong* PICC) or may be situated in the **external** end of the catheter (eg the *PASV* PICC).
- **Patients may return home with a PICC in situ**, and therefore patient education regarding the recognition and reporting of complications is of great importance. The PICC usually exits onto the patient's arm and so it is generally not practical for the patient to care for the catheter him/herself. Liaison with the primary health-care team is vital.
- **Placement is contraindicated** following axillary node dissection or irradiation, or in the case of lymphoedema of the arm, axillary node disease or skin infection at the insertion site¹⁵.
- **The PICC should not be confused with a "midline catheter"** which is usually "20cm in length, with the tip terminating in the region of the axillary vein, and is designed for short-term peripheral drug delivery"¹⁶. A midline catheter is not a Central Venous Catheter.

d) Implantable Ports *Often called Portacaths*

- See Chapter 4 Care of Implantable Ports for care of these catheters.
- **The Implantable Port is similar to a tunnelled line** but instead of protruding from the patient's chest, the catheter terminates in a self-sealing injection port which is implanted under the skin. **There are therefore no external parts.** The port is accessed through the skin using a dedicated non-coring needle (see diagram next page)
- Some patients find an Implantable Port **more discreet and less intrusive** than a tunnelled CVC¹⁷. **Ports require less maintenance** when not in use than other types of catheter. They may also offer a lower risk of infection^{17,18}.
- Implantable Ports are suitable for patients who require **long-term frequent and intermittent venous access**. Arguably they are less than ideal for long-running continuous infusions because of the risk of needle dislodgement¹⁹. **The patient may return home with the port in situ**, and therefore patient education regarding the recognition and reporting of complications is of great importance.

Diagram 6. Implantable Port

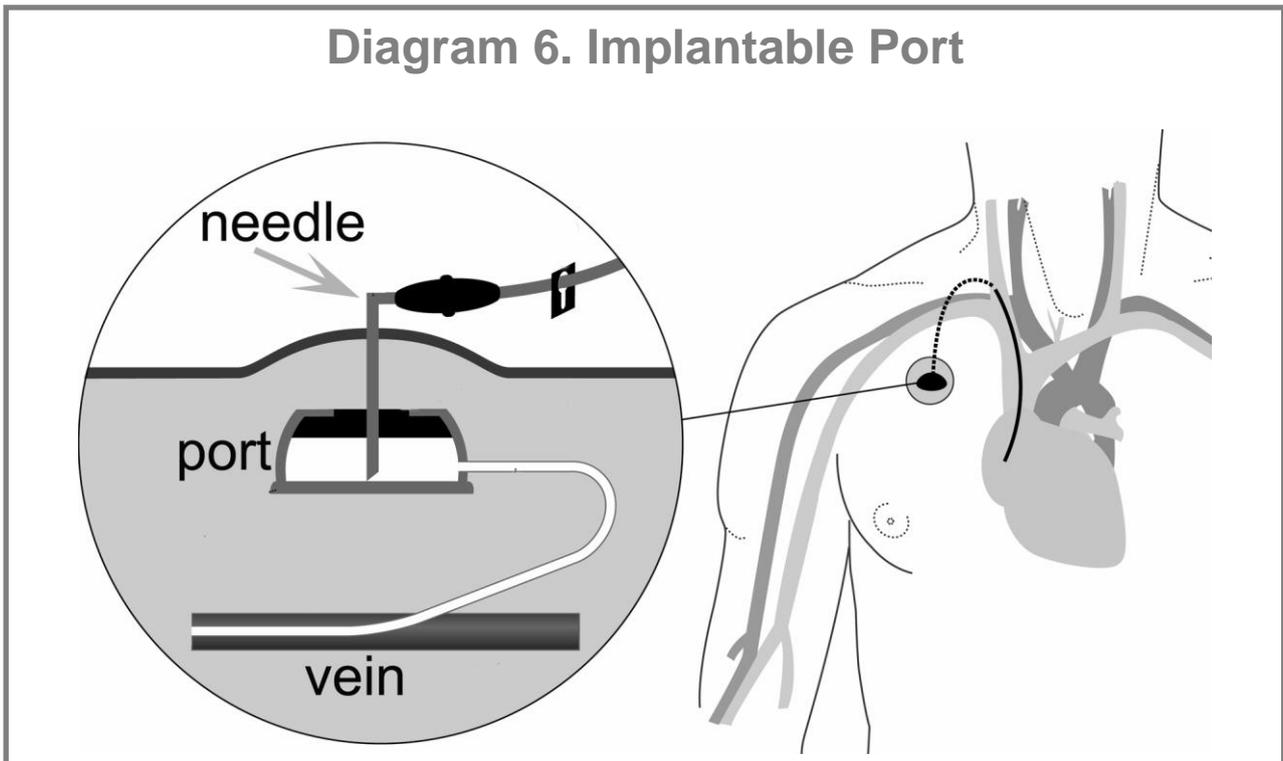

- **Dual lumen devices are available.** These are equipped with two access ports side-by-side which can be accessed separately using two different needles. Each lumen provides independent access to the venous circulation, so that incompatible drugs/fluids may be administered simultaneously.
- Ports may also be used as **an alternative to subcutaneous administration** of long-term maintenance therapies when the subcutaneous route has become unacceptable to the patient or unreliable – e.g. due to subcutaneous nodule formation.
- **Placement is not recommended** in obese or cachexic patients, before or after chest irradiation, or at mastectomy sites¹⁷.

e) CVCs used for Blood Processing (eg Apheresis) Often called *Permacaths* / *Vascaths*.

- See Chapter 5 Care of CVCs Used For Blood Processing for care of these catheters.
- CVCs used for blood processing – eg apheresis - are very similar to the catheters described in a) and b) above. **They can be non-tunnelled (eg *Vascaths*) or tunnelled (eg *Permacaths*).**
- **These catheters differ from other CVCs in the following respects:**
 - **Larger lumen size** compared to other CVCs.
 - **The internal tip of the catheter is designed differently** so as to allow blood to be withdrawn freely via one lumen and returned via the other lumen downstream of the blood being withdrawn (thus avoiding recycling of the treated blood). Confusingly, the lumens are often colour-coded red and blue and referred to as the “arterial” and “venous” lumens. In fact both lumens lead into a vein and not an artery.

- **In some settings these catheters are locked between uses with an exact volume of concentrated heparin** solution to minimise the risk of occlusion. This varies depending on the patient's clinical status and local guidelines. If heparin is used as a lock, the lumen must be clearly labelled and the heparin must be withdrawn from the catheter before use otherwise the patient will receive an unacceptably high dose of heparin. See Chapter 5 Care of CVCs Used for Blood Processing for details.
- **In other respects these catheters are identical** to the CVCs described in a) and b) above.

Appendix II: Background Information: Principles of Care

(i) General Principles

- A recommendation of the EPIC-3 guidelines is that healthcare workers should be aware of the manufacturer's advice relating to individual catheters, connection and administration set dwell time, and compatibility with antiseptics and other fluids to ensure the safe use of devices.
- **Infection** is the most common complication associated with central venous access and one of the most serious. You must be familiar with infection control guidance when caring for patients with CVCs. Refer to your local hospital's policy.
- **Decontaminate hands** before and after each patient contact using correct hand hygiene procedure.
- **Use an aseptic, non-touch technique (ANTT) following** whenever the CVC is accessed and during procedures involving exit sites. *To prevent infection. A strong correlation exists between bacteraemia and the presence of a CVC²⁰.*
- **Regularly inspect for signs of infection** (at least daily if patient is in hospital) including inspection of exit site and monitoring of temperature, pulse and blood pressure at least daily when the patient is in hospital. *To detect infection*
- **Take action immediately** if there are signs of CVC-related infection. These include:
 - Pyrexia, rigor, malaise (see Chapter 6 Managing Complications)
 - Tenderness, inflammation and or pain at exit site (see Chapter 6 Managing Complications)
- **Wear gloves when carrying out dressing changes and when accessing the catheter.** Whether or not these should be sterile will depend on the task. *Gloves should be worn to prevent descaling of bacteria onto key parts²¹.*
- **The practice of administering prophylactic antibiotics** at the time of CVC insertion should NOT be routinely followed. (The exception is if the patient's MRSA status is unknown or positive, in which case follow current UCLH guidelines or your local trusts guidelines in this circumstance). *The Department of Health's Epic Guidelines on the prevention of infection in Central Venous Catheters (phase 2) specifically states that this practice is not supported by research and may encourage resistant organisms⁹.*
- **Do not allow air to enter the catheter.** The catheter should never be left open to air entry and all syringes and intravenous administration sets must be carefully primed. *To prevent air embolism. The negative pressure within the chest may suck air into the catheter during inspiration especially if the patient is sitting up²².*
- **Cap off the catheter with a needle-free access device** (eg *Bionector*[®]) when not in use. This will minimise interruptions to the closed system. Unless manufacturer's instructions vary, this should be changed every 7 days or every 100 uses, whichever is the sooner. *Risk of contamination increases with every interruption to the closed system²⁰.*
- **When accessing the catheter through the needle-free device, always decontaminate the device using Chlorhexidine 2% in alcohol and using ANTT (see (ii) Accessing the Catheter below).** The device should be cleaned each time using a 30-second vigorous friction rub and allowed to dry before inserting a sterile syringe-tip.

- **Whenever the bung/access device is removed** from the catheter then it must be replaced with a new, sterile device. *To prevent infection.*
- **If the catheter possesses an integral clamp, keep it closed** whenever the cap is removed and at all other times except when administering or withdrawing fluids. Clamping should always take place at the designated area and never at the thickened area near the hub. *The clamp will prevent air entry and bleeding should the luer lock cap become unattached. Repeated clamping away from the specially reinforced area may result in damage to the catheter.*
- **The practice of administering prophylactic mini-dose Warfarin** to patients with CVCs should NOT be followed. *Mini-dose Warfarin has recently been shown to be ineffective in the prevention of thrombosis in cancer patients with CVCs²³. (NB dose adjusted Warfarin did show some efficacy but with an increased risk of serious bleeding).*
- **Should the catheter fracture or be accidentally cut**, clamp it atraumatically without delay proximal to the break. Specialist advice should be sought immediately to consider removal or repair of the catheter. *To prevent haemorrhage, air embolism and infection.*
- **Always secure the catheter firmly to the skin.** *For patient's comfort, to prevent tension or accidental dislodgement, and to reduce 'to and fro' motion which increases the risk of catheter related sepsis²⁰.*
- **Central Venous Catheters** should be removed as soon as possible if they are not needed.

(ii) Accessing the Catheter

- **Decontaminate hub using Chlorhexidine 2% in alcohol (eg Clinell wipe) using a 30-second friction rub and allow to dry before accessing any lumen. Use an aseptic, non-touch technique (ANTT) following local guidelines** whenever the CVC is accessed. *To prevent bacteria from entering the blood stream. A strong correlation exists between bacteraemia and the presence of a CVC²⁴.*
- **Before it is used for administering therapeutic drugs or fluids, the patency and correct functioning of the catheter should be established²⁵.** Signs of catheter occlusion, whether partial or complete, should be taken seriously and action should be taken earlier rather than later to restore full patency. Ignoring the early signs may lead to the development of more serious problems which cannot then be easily rectified – eg complete blockage²⁶. In addition unresolved patency problems may mask a malpositioned catheter which carries and increased risk of thrombosis.
- **Nurses using CVCs can be confident of access** if all three of the following apply
 - The catheter can be **flushed with ease**.
 - **Blood can be withdrawn** from the catheter (see below for technique)
 - The patient experiences **no discomfort** during flushing/infusion and there are no other complications
- **If any of these criteria are not met** see Chapter 6 Managing Complications.
- **Ways of assessing these three criteria will vary with the setting.** Here are some points to note:
 - **A proper assessment of the catheter involves observing the exit site and the area around** as this may reveal any signs of thrombosis, leakage, infection etc. While this is not

necessarily appropriate every time the catheter is used it should be a regular part of your practice.

- **When checking for flashback of blood, do not discard blood unnecessarily. To assess for flashback you can either**
 - attach a syringe containing 10mls 0.9% saline (must be prescribed) to the catheter, flush a couple of mls into the line and then withdraw. As soon as you see a trace of blood in the catheter or syringe just flush the rest of the saline into the line using the push-pause technique as described in (iii) a) below.
 - or use a gravity technique (ie with clamps open briefly hold an attached infusion below the level of the patient's heart until you see flashback of blood). NB this may not work for valved PICCs (eg PASV or Groshong).
 - BUT note that if there are infusional vasoactive drugs in the lumen, withdraw prior to flushing to avoid bolus dose.
- **For patients receiving parenteral nutrition (PN)**, testing for flashback via the PN lumen on a regular basis is not advised as it may increase the infection risk¹². But note that a full assessment of any central venous catheter is advised if there are doubts about patency or if malposition is suspected.
- **Assessing CVCs in patients requiring blood processing (eg apheresis) requires specialist knowledge:** refer to Chapter 5 Care of CVCs used for Blood Processing for care of these patients.

(iii) Flushing After and Between Uses

(a) Flushing Technique:

- **Where possible, do not use syringes smaller than 10 mls** for infusion into the catheter¹⁵. *To prevent excessive pressure being exerted on the lumen which might cause it to rupture. Smaller syringes exert greater pressure. But please note that syringe size alone is not sufficient to prevent rupture. "When resistance is felt, if more pressure applied to overcome it, catheter fracture could result regardless of the syringe size..."*²⁶.
- **Use a brisk 'push-pause' flushing technique** routinely when flushing the catheter - i.e. flush briskly, pausing briefly after approximately each ml of fluid. *The 'push-pause' technique causes turbulence within the catheter, which helps to flush away any debris and prevent occlusion of the lumen*^{15,25}.
- **If the catheter possesses a clamp, clamp the line while the final ml of the flush is being injected.** If there is no clamp you can achieve a "positive pressure finish" by removing the syringe from the *Bionector*[®] (or similar) while injecting the last ml, being careful to avoid any spray from the syringe. *Maintaining positive pressure helps prevent blood entering the catheter after flushing, which might lead to occlusion or thrombus formation*¹⁵.
- **Do not routinely withdraw and discard blood from the catheter before flushing** in an attempt to avoid flushing bacteria and clots into the patient²⁵. *There is no evidence that withdrawing prior to flushing reduces infection or embolism. But note that if the catheter is to be used for administering drugs or fluids, checking for "flashback" should be a routine part of catheter assessment:* see ii) Accessing the Catheter (above).

- **For patients receiving parenteral nutrition (PN)**, testing for flashback via the PN lumen on a regular basis is not advised as it may increase the infection risk¹². But note that a full assessment of any central venous catheter is advised if there are doubts about patency or if malposition is suspected.

(b) Frequency of flushing and flushing solutions:

- This varies depending on the device. See Chapters 1 – 5.
- *Hepsal* (heparinised saline 10 u/ml) should no longer be used to flush in PICCs and Hickman lines. This is due to the lack of convincing evidence to support its use coupled with the increasing concerns raised at a national level about the potential dangers of heparin.
- Instead of the familiar “10mls 0.9% saline plus 5mls Hepsal” combination for locking off these lines, we are now advocating 10mls 0.9% saline *and repeat*. The rationale behind repeating the 10mls 0.9% saline is that using just 10mls instead of the 15mls involved in the saline/Hepsal combination may lead to more occlusion problems. (Indeed this was anecdotally felt to be case when one ward changed their practice this year.) The rationale for using two 10ml syringes rather than a 20ml syringe is that a 10ml syringe creates higher pressures and therefore greater turbulence in the line (which is important in maintaining line patency).
- But note that heparin / heparinised saline continues to be used in implantable ports and CVCs used for blood processing.

(iv) Care Of The Exit Site

a) Dressings:

- **Patients aged 16+:** *Tegaderm™ CHG* or equivalent dressing changed every 7 days or sooner if dressing becomes detached or soiled. *Tegaderm™ CHG has been shown to be highly effective at reducing skin flora for up to 10 days*²⁷. A recommendation of the EPIC-3 guidelines is to use a chlorhexidine impregnated sponge dressing at the exit site as a strategy to reduce catheter related blood stream infection.⁵⁵
- **Patients under 16: IV-dedicated occlusive transparent dressing**, changed every 7 days²⁵. *Some researchers have found iv-dedicated transparent dressings to be associated with a lower risk of infection than other dressings*²⁸.
- **Patients who cannot tolerate *Tegaderm™ CHG*:** use a plain iv dedicated transparent dressing (plus *Biopatch®* for patients aged 16+) changed every 7 days.
- **Patients who cannot tolerate a transparent dressing at all:** use a sterile dry dressing (plus *Biopatch®* for patients aged 16+), changed at least weekly or more frequently if necessary for inspection or adherence²⁵.
- **No dressing.** This may be suitable for some patients with tunnelled CVCs from 21 days post insertion once the tissues have fibrosed around the cuff and in the absence of exudate or signs of infection. *“No dressing” performed just as well as 3 types of dressing in one study comparing infection rates*^{20, 29, 30}.
- As a general principle, where a dressing is used it should be **inspected regularly and renewed immediately should it become soiled, wet or detached**²². *A moist environment is one in which bacteria readily multiply*³¹.

- **If the exit site is reddened, painful, exudating or infected**, take a swab and inform the medical team with a view to antibiotics and / or removal of the device. Stop using *Tegaderm™ CHG* and *Biopatch®*. Increase the frequency of dressing change depending on the amount of exudate.

b) Dressings Immediately post insertion:

- Exit sites may bleed immediately post insertion so it may not be appropriate to apply *Tegaderm™ CHG* or **equivalent** until bleeding has stopped. In this case a gauze pad covered with an iv-dedicated transparent may be used.

c) Cleaning of Exit Site:

- **Cleaning should be carried out using ANTT.**
- **Clean exit site at dressing changes** with a sterile single-use swabstick containing chlorhexidine 2% in alcohol (eg *Chloraprep Sepp®* or *Chloraprep 3mls®*) and using a 30-second back and forth friction rub. Allow to dry. (NB *Chloraprep®* can be used for patients aged 2 months and over.)
- **Loose blood, exudate or other debris** which might provide a focus of infection or might impair inspection of the wound may be gently removed with sterile gauze soaked with 0.9% saline* prior to cleaning with Chlorhexidine³².

(v) Removal

Some CVCS are simple and relatively safe to remove. With others, there is high risk of air embolism³³ and so removal requires a higher level of training and skill. See Chapters 1 – 5 for guidelines on removal.

Appendix III: Background Information: Complications of CVCs

(i) Pneumothorax

A pneumothorax is the presence of air in the pleural space between the lungs and the chest wall. It can occur during the insertion of a CVC when a needle used to access the subclavian or jugular veins inadvertently punctures the lung. The risk of pneumothorax is considerably reduced when ultrasound guidance is used to access the veins. Since most CVCs used in Cancer Services are inserted in this way, a chest x-ray to screen for pneumothorax following insertion is no longer routinely carried out. If the catheter has been inserted using fluoroscopy to screen the tip position, then there will be no need for a post insertion chest x-ray at all. However if the catheter has been inserted without fluoroscopy a post-insertion chest x-ray will still be required to check the tip position. The exception is with short term femoral CVCs. Here there is no risk of pneumothorax and the tip position is not routinely checked.

A pneumothorax may be clinically silent, or may lead to a life-threatening emergency situation with respiratory distress, reduced oxygen saturation levels, tachycardia or hypotension. A small pneumothorax may resolve spontaneously. In severe cases a chest drain may be necessary.

(ii) Infection

Infection is the most common complication associated with central venous access³⁴ and one of the most serious³⁵ with estimated mortality rates ranging from 1 – 35%³⁶.

Refer to local care guidelines e.g. Central Venous Catheter Care Bundle to document care around these lines.

Contamination can occur during insertion of the CVC or at a later stage via the hands of healthcare workers, or transferred from the patient's skin or other anatomical sites. Infection may be relatively minor or may be life-threatening.

Bacteria can colonise a CVC either on its exterior or interior surface: ie colonisation is either **extraluminal** or **intraluminal**. **Extraluminal** infections usually begin at the exit site and may remain confined to that area or may track along the catheter into the bloodstream. **Intraluminal** infections are caused by contamination via the hub of the catheter³⁷.

A Microbiology opinion should usually be sought in deciding the best action to take in the event of signs of infection. Exit site infections can often be treated successfully with antibiotics, especially in PICCs, and in tunnelled CVCs (Hickman lines) where the vein and the exit site are separated by the tunnel. In non-tunnelled centrally inserted CVCs, however, treatment is less likely to be successful, as there is less distance between the exit site and the blood stream⁹. By the same token, infections in tunnelled CVCs involving the skin tunnel itself above the cuff are notoriously difficult to treat and the same applies in implantable ports where there is infection of the port pocket.

The risk of infection can be reduced by strict adherence to best practice e.g. UCLH Central Venous Catheter Bundle and using Aseptic Non-Touch Technique. Intravenous tubing and stopcocks should be changed according the UCLH Infection Control Guidelines or according to your local policy . If Parenteral Nutrition is to be given, one lumen should be used exclusively for this purpose.

(iii) Thrombosis

Thrombosis occurs when a clot develops within the vein around the catheter. Unless the clot is at the internal tip of the catheter, it will not usually affect the patency of the catheter. Thrombosis formation is a natural response to vascular injury. Damage to the vessel wall can occur during catheter insertion, or may be due to mechanical or chemical irritation in an incorrectly placed catheter eg where the tip of the catheter is in too small a vein, or rubbing against the vein wall instead of floating parallel to it.

The risk of thrombosis is increased in patients who are pregnant or immobile or who have diabetes or cancer. Surgery, chemotherapy, hormonal agents, haemodialysis and CVC-related infection are all thought to be risk factors^{11,26,38,39}. It used to be thought that minidose Warfarin might reduce the risk of thrombosis in Cancer patients, but this has recently been disproved²³.

Patients who develop thrombosis are at increased risk of pulmonary embolism and infection³⁸.

A large proportion of cancer patients with CVCs have thromboses which are never detected^{40,41}. When a thrombosis does become symptomatic, it will usually cause swelling of the arm, neck, face or lower limb. There may be associated pain, tingling or numbness, distended neck or peripheral veins^{42,43}. The presence of a thrombosis can usually be confirmed by use of Doppler ultrasound.

Clinically evident thrombosis is more common in cancer patients with PICCs than with tunnelled CVCs or implantable ports, probably because they occupy smaller veins which are more easily occluded.

Unless the catheter is incorrectly positioned, it may be possible to treat a thrombosis using anticoagulants without removing the catheter.

(iv) Air Embolism

An air embolism is a potentially fatal complication. It can happen at any stage if air is allowed to enter the catheter – eg if a catheter is left unclamped when the cap is removed – but is most likely to occur during the insertion or removal of the catheter. The risk is increased if the patient is dehydrated, is unable to lie flat, or has an uncontrolled cough at the time of insertion or removal.

As with pneumothorax, air embolism may be clinically silent or may be accompanied by any or all of the following: anxiety, cyanosis, dyspnoea, tachycardia, hypotension, chest pain, loss of consciousness and death.

(v) Cardiac Arrhythmias

Atrial or ventricular arrhythmias can occur when the tip of the CVC is placed within the heart^{44,45}. In practice, non-tunnelled and tunnelled CVC tips correctly placed in the right atrium rarely cause arrhythmias. PICCs are probably most likely to cause problems because the PICCs are more floppy and more likely to move within the vein than other CVCs. They can also move further into the heart as the patient moves his / her arm. Arrhythmias caused in this way will usually resolve when the catheter is pulled back by a few centimetres. Any patient experiencing unresolved palpitations or arrhythmias should be assessed by a medical team as soon as possible.

(vi) Cardiac Tamponade

This is a rare complication of CVCs, seen mainly in neonates. Cardiac tamponade arises when fluid (in this case blood) accumulates in the pericardial space around the heart and impairs cardiac function. This is a

catastrophic, often fatal event. The patient is likely to exhibit a sudden onset of severe cardiorespiratory symptoms.

Cardiac tamponade can arise in a patient with a CVC if the heart is punctured either during insertion or subsequently by a malpositioned catheter.

(vii) Patency Impairment

Patency is said to be impaired in any of the following situations:

- The catheter is completely blocked and cannot be flushed at all.
- The catheter can be flushed using a syringe but there is sluggish, absent or intermittent free-flow when infusion of fluids by gravity is attempted.
- The catheter flushes easily but aspiration of blood is sluggish or absent.

Patency problems should be taken seriously. Ignoring the early signs may lead to the development of more serious problems which cannot then be easily rectified – eg complete blockage or thrombosis²⁶.

The causes of patency problems include

Clotted blood within the catheter: This can be avoided by good flushing techniques as described in these guidelines. When problems do arise, they can usually be solved relatively easily by use of a thrombolytic such as Urokinase: see Chapter 7 Maintaining Patency.

Fibrin Sheath: Fibrin sheaths are thought to occur in most CVCs left in place for over 7 days⁴⁶. A fibrin sheath is a kind of sleeve made of a fibrous collagen substance which can form around the catheter within the blood stream. It may extend to form a kind of “sock” protruding beyond the tip of the catheter, and if this happens it may impair the patency of the catheter: most commonly it will prevent blood from being withdrawn from the catheter because the fibrin sheath is sucked against the tip of the catheter. In severe cases a fibrin sheath may also lead to backtracking of infused fluids between the fibrin sheath and the catheter, causing leakage of those fluids into the tissues⁴⁷. Fibrin sheaths may be diagnosed using fluoroscopy and are associated with an increased risk of infection⁴⁸ as they provide an ideal medium for the proliferation of bacteria. They can sometimes be removed by stripping under imaging or by an infusion of a thrombolytic.

Mechanical obstruction: A mechanical obstruction can occur internally or externally. Internal obstruction may be due the catheter being incorrectly positioned: eg it may be kinked or the tip of the catheter may be resting against a vessel wall rather than floating free within the bloodstream (see Incorrect Position below). This might be due to poor insertion technique, or due to catheter dislodgement. A simple chest x-ray will often reveal an incorrectly positioned catheter. External kinking of the catheter can also cause patency problems: its’ worth checking for a bra-strap or an over-tight stitch before looking for a more complicated cause!

Build up of lipids from parenteral nutrition or drug precipitation within the catheter caused by too high a concentration or incompatibility of drugs: If this appears to be a likely cause of occlusion, consult Medical/Pharmacy advice for a suitable agent to dissolve occlusion.

(viii) Incorrect Position

A CVC should be considered to be in an incorrect position when any of the following apply:

- The tip is not in the right atrium, the superior or the inferior vena cava.

- The tip of the catheter is not floating freely parallel to the vein wall.
- The catheter is kinked within the body or pinched between internal structures

Incorrect position may be the result of poor insertion technique or may occur spontaneously in a previously well-positioned catheter. It is not unknown for a CVC to “migrate” within the venous system for no apparent reason. Hadaway reports that “Changes in intrathoracic pressure, coughing, sneezing, Valsalva manoeuvre such as during heavy lifting, vigorous extremity use, forceful flushing, or congestive heart failure could lead to migration of the tip”²⁶. In addition the catheter may become dislodged if it is not correctly secured in place, or is accidentally pulled.

If a CVC is incorrectly positioned there is a high risk of thrombosis and patency impairment⁴⁹. If it is kinked internally there is also the risk that the catheter may split, leading to extravasation of drugs / fluids and in serious cases, embolisation of the catheter itself.

You should suspect incorrect position if there are patency problems despite the use of a thrombolytic, if the patient complains of pain on flushing, if the external length of the catheter increases, if the patient develops a thrombosis, or if the cuff of a tunnelled CVC protrudes from the exit site⁵⁰.

A malpositioned, kinked or pinched catheter should be repositioned, replaced or removed as soon as practicable (except PICCs in certain situations discussed below). Leaving it in place for any length of time represents a high risk of thrombosis and/or catheter fracture / embolism.

Immediately following insertion, PICCs are sometimes found on X-ray to have fed up into the jugular vein, across into the opposite subclavian, or back down an arm vein. In these cases it may be worth leaving the PICC in overnight or flushing briskly with 20ml 0.9% saline* and then repeating the X-ray as the PICC will often move into the Superior Vena Cava^{51, 52}. Discuss with the person inserting the PICC and patient’s medical team.

(ix) Extravasation of Fluids / Drugs due to Incorrect Needle Position or Needle Dislodgement (in Implantable Ports)

The non-coring needle should be correctly placed into the port (see Appendix I (iv) d) Diagram 6). If the needle is not inserted far enough into the port or if the needle misses the port altogether this may lead to fluids / drugs being infused into the subcutaneous tissues.

The needle may become dislodged if it is inadequately secured with dressing tape, if there is tension on the extension tubing or if the needle used is of insufficient length, causing the patient’s normal movements to loosen the needle. The problem will usually be noticed when there is discomfort and/or oedema at the entry site combined with lack of free-flow of fluids.

If extravasation has occurred or is suspected, the needle should be removed and a fresh needle used to access the port correctly. If vesicant or irritant solutions (eg chemotherapy) are extravasated, seek medical / pharmacy advice and refer to the local trust Cytotoxic Policy⁵³.

(x) Catheter Fracture

This may occur **externally** or **internally** and may result from over-forceful flushing, trauma to the catheter or incorrect position (eg kinking leading to wear-and-tear).

An **external** fracture will result in leakage of blood or fluids from the catheter. Sometimes there is an obvious fracture. The line must be clamped or folded over on itself immediately to prevent air embolism. Sometimes the catheter can be repaired or replaced over a guidewire but the advisability of this will depend on the patient's clinical status. In addition, unless the correct equipment and expertise are available for a repair to be carried out, the catheter should be removed immediately, as there is a high risk of infection and air embolism.

Internal fracture will usually result in patency impairment and / or pain, redness and swelling when the catheter is flushed.. There is a risk that the catheter itself will embolise. If this occurs there may be no symptoms at all or there may be signs of pulmonary embolism. ie acute onset of any or all of the following - anxiety, pallor, cyanosis, shortness of breath, rapid weak pulse, hypotension, chest pain, loss of consciousness.

(xi) Separation of port and catheter (in Implantable Ports)

This is rare⁵⁴ but should always be considered when problems arise with patency of the port or there is Extravasation with associated discomfort and oedema despite proper position of needle.

As with catheter fracture (see (x) above) there is a risk that the catheter may embolise. Surgical removal or repair of the port and catheter is essential if separation is confirmed.

References

- ¹ Mayo DJ "Administering Urokinase: Clearing the Way" Nursing98 December
- ² Vesely, T. "Central Venous Catheter Tip Position: A Continuing Controversy" Journal of Vascular and Interventional Radiology, Volume 14(5) May 2003, pp 527-534
- ³ Racadio, JM, Doellman DA, Johnson ND, Bean JA, Jacobs BR. "Pediatric Peripherally Inserted Central Catheters: Complication Rates Related To Catheter Tip Location." Pediatrics. 107(2):E28, 2001 Feb
- ⁴ Puel V et al 1993. "Superior vena cava thrombosis related to catheter malposition in cancer chemotherapy given through implanted ports". Cancer. 72(7):2248-52, 1993 Oct 1
- ⁵ Eastridge BJ and Lefor, AT. "Complications of indwelling venous access devices in cancer patients". Journal of Clinical Oncology. 13(1):233-8, 1995 Jan.
- ⁶ Schuster M, et al. "The carina as a landmark in central venous catheter placement". British Journal of Anaesthesia 2000; 85: 192–4
- ⁷ Fletcher S, Bodenham AR. "Safe placement of central venous catheters: where should the tip of the catheter lie?" British Journal of Anaesthesia 2000; 85: 188–91.
- ⁸ Department of Health Saving Lives Campaign: <http://www.clean-safe-care.nhs.uk>
- ⁹ Pratt RJ et al (2007) "epic2: National Evidence-based Guidelines for Preventing Healthcare-Associated Infections in NHS Hospitals in England". Journal of Hospital Infection 2007; 65(Supplement):S1-S64
- ¹⁰ National Institute for Clinical Excellence (September 2002) "Guidance on the Use of Ultrasound Locating Devices for Placing Central Venous Catheters." NICE Technology Appraisal No 49. London: National Institute for Clinical Excellence. Available from www.nice.org.uk
- ¹¹ Grove, Jay R and Pevec, William C "Venous Thrombosis Related to Peripherally Inserted Central Catheters" Journal of Vascular and Interventional Radiology Volume 11(7) July/August 2000 pp 837-840
- ¹² UCL Hosp NHS Trust "Adult and Adolescent Parenteral Nutrition Handbook" 2008 <http://insight/guidelines/Trustwide%20Guidelines/Adult%20and%20Adolescent%20Parenteral%20Nutrition%20Hand-book.doc>
- ¹³ Wickham R et al "Long-term CVCs - Issues for Care" Seminars in Oncology Nursing Vol 8 No 2 May 1992 pp 133-147
- ¹⁴ Wilson J A "Preventing Infection During IV Therapy" Professional Nurse March 1994 pp 338-392
- ¹⁵ Todd J "Peripherally inserted central catheters" Professional Nurse Vol 13 No 5 Feb 1998 pp 297-302
- ¹⁶ Todd, J "Peripherally Inserted Central Catheters and their Use in IV Therapy" British Journal of Nursing Vol 8 No 3 1999 pp 140-48

- ¹⁷ Camp-Sorrell D "Implantable Ports - Everything you Always Wanted to Know" Journal of Intravenous Nursing Vol 15 No 5 Sep/Oct 1992 pp 262 - 272
- ¹⁸ Krzywda, E et al "Catheter Infections: Diagnosis, Etiology, Treatment, and Prevention" Nutrition in Clinical Practice Vol 14 No 4 August 1999 pp178 - 90
- ¹⁹ Schulmeister L "Needle Dislodgement from Implanted Venous Access Devices; Inpatient and Outpatient Experiences" Journal of Intravenous Nursing Vol 12 No 2 March/April 1989 pp 90-92
- ²⁰ Haller L and Rush K "CVC infection: a review" Journal of Clinical Nursing Vol 1 1992 pp 61-66
- ²¹ Rowley, S "Aseptic Non Touch Technique (ANTT)" Nursing Times Feb 15th Vol 97 No 7 2001: Infection Control Supplement V1-V111
- ²² Cornock M "Making Sense of CVCs" Nursing Times Vol 92 No 49 Dec 4th 1996 pp 30-31
- ²³ Young, A "WARP - A multicentre prospective randomised controlled trial (RCT) of thrombosis prophylaxis with warfarin in cancer patients with central venous catheters (CVCs)" 2005 ASCO Annual Meeting
- ²⁴ Haller L and Rush K "CVC infection: a review" Journal of Clinical Nursing Vol 1 1992 pp 61-66
- ²⁵ RCN IV Therapy Forum "Standards for Infusion Therapy" Royal College of Nursing July 2003
- ²⁶ Hadaway L "Major Thrombotic and Nonthrombotic Complications: Loss of Patency" Journal of Intravenous Nursing Vol 21 No 5S September/October 1998
- ²⁷ Maki DG et al; "A Novel Integrated Chlorhexidine-impregnated Transparent Dressing for Prevention of Vascular Catheter-related Infection: A Prospective Comparative Study in Healthy Volunteers". SHEA Annual Scientific Meeting 2008 (Poster).
- ²⁸ Maki D G et al "Prospective Randomised Trial of Povidone Iodine, Alcohol and Chlorhexidine for Prevention of infection Associated with Central Venous and Arterial Catheters" Lancet 383 1991 pp 339-343
- ²⁹ Krzywda, E "Predisposing Factors, Prevention and Management of Central Venous Catheter Occlusions" Journal of Intravenous Nursing Vol 22 No 6S November/December 1999 pp S11 – S17
- ³⁰ Olson, K et al "Evaluation of a No-dressing Intervention for Tunnelled Central Venous Catheter Exit Sites" Journal Of Infusion Nursing Volume 27(1) January/February 2004 pp 37-44
- ³¹ Gabriel J "Care and management of peripherally inserted central catheters" British Journal of Nursing Vol 5 No 10 1996 pp 594-599
- ³² Oliver L "Wound Cleansing" Nursing Standard Vol 11 No 20 Feb 5th 1997 pp 47-51
- ³³ Drewett, S "Central venous catheter removal: Procedures and rationale" British Journal of Nursing. London: Dec 8, 2000-Jan 10, 2001.Vol.9, Iss. 22; pg. 2304

- ³⁴ Krzywda, Elizabeth A “*Central Venous Catheter Infections: Clinical Aspects of Microbial Etiology and Pathogenesis*”. Journal of Infusion Nursing Volume 25(1) January/February 2002 pp 29-35
- ³⁵ Oncu, S and Sakarya, S. “*Central venous catheter-related infections: an overview with special emphasis on diagnosis, prevention and management*”. Internet Journal of Anesthesiology, 2003, vol. 7, no. 1 .
- ³⁶ Citton-R et al. “*Old and new tools in the diagnosis of central venous catheter-related bloodstream infections: Is there a role for brushing?*” Journal of Vascular Access 2004 Volume 5 Issue 1 Pg 10-12
- ³⁷ Hall, K and Farr, B. “*Diagnosis and Management of Long-term Central Venous Catheter Infections*”. J Vasc Interv Radiol 2004; 15:327–334
- ³⁸ Marinella, Mark A et al “*Spectrum of upper-extremity deep venous thrombosis in a community teaching hospital*” Heart and Lung: The Journal of Acute and Critical Care Volume 29(2) March/April 2000 pp 113-117
- ³⁹ Lee, Agnes Y and Levine, Mark N. “*Management of Venous Thromboembolism in Cancer Patients*” Vol 14, no 3, (March 2000)
- ⁴⁰ Balestreri-L et al “*Central venous catheter-related thrombosis in clinically asymptomatic oncologic patients: A phlebographic study.*” European Journal of Radiology {EUR-J-RADIOL}, 1995, Vol/Iss/Pg. 20/2 (108-111).
- ⁴¹ Allen, Anthony W et al. “*Venous Thrombosis Associated with the Placement of Peripherally Inserted Central Catheters*”. Journal of Vascular and Interventional Radiology Volume 11(10) November/December 2000 pp 1309-1314
- ⁴² Jacobs P et al “*Chest Pain As The Presenting Symptom In Catheter-Associated Thrombosis Of The Superior Vena-Cava*” S Afr Medical Journal 88(10) 1998 pp 1284-5
- ⁴³ Kayley, J “*Skin-Tunnelled Cuffed Catheters*” Community Nurse June 1997 pp 21-22
- ⁴⁴ Bivins MH and Callahan MJ “*Position-Dependent Ventricular Tachycardia Related To A Peripherally Inserted Central Catheter*” Mayo Clinic Proceedings, 75 (4): 414-6, 2000 Apr
- ⁴⁵ Teichgräber, UK et al “*Central Venous Access Catheters: Radiological Management of Complications*” Cardiovascular and Interventional Radiology (2003) 26:321-333
- ⁴⁶ Gabriel J “*Fibrin sheaths in vascular access devices*” Nursing Times Vol 93 No 10 March 5 1997
- ⁴⁷ Mayo DJ and Pearson DC “*Chemotherapy Extravasation: A Consequence of Fibrin Sheath Formation Around Venous Access Devices*” Oncology Nurse Forum Volume 22 No 4 May 1995, 675-680
- ⁴⁸ Mehall, JR, Saltzman DA, Jackson RJ and Smith SD “*Fibrin Sheath Enhances Central Venous Catheter Infection*” Critical Care Medicine Volume 30 (4) April 2002, 908-912

⁴⁹ Aitkin DR and Minton JP "The 'Pinch-off Sign;' A Warning of Impending Problems with Permanent Subclavian Catheters" American Journal of Surgery Vol 148 Nov 1984 pp 633-636

⁵⁰ Jones GR "A Practical Guide to Evaluation and Treatment of Infections in Patients with Central Venous Catheters" Journal of Intravenous Nursing Vol21 No 5S September/October 1998 pp S134 – S 142

⁵¹ Banks N "*Positive Outcome after Looped Peripherally Inserted Central Catheter Malposition*" Journal of Intravenous Nursing Vol 22 No 1 January/February 1999 pp 14 - 18

⁵² Rastogi S et al "Spontaneous Correction Of The Malpositioned Percutaneous Central Venous Line In Infants" Pediatric Radiology. 28(9): 694-6, 1998 Sep

⁵³ North London Cancer Network / North East London Cancer Network "*Guidelines for the Safe Prescribing, Handling and Administration of Cytotoxic Drugs.*" Latest version at time of writing: January 2005

⁵⁴ Moore C et al "*Nursing Care and Management of Venous Access Ports*" Oncology Nursing Forum Vol 13 No 3 May/June 1986 pp 35-39

⁵⁵ H.P. Loveday et al epic3: National Evidence-Based Guidelines for Preventing Healthcare-Associated Infections in NHS Hospitals in England Journal of Hospital Infection 86S1 2014 S1–S70